



PRODUCT CATALOGUE

The choice of dental professionals worldwide.

TABLE OF CONTENTS

CLINICAL RA - STEEL SHANK MULTI-USE

PAGE	PRODUCT	APPLICATION	FLEXIBILITY
MULTIPLE MATERIALS			
8	KENDA NOBILIS	For Composites, Ceramics & Hybrid-Ceramics	Flexible
9	KENDA UNICUS	For Composites, Ceramics & Hybrid-Ceramics	Flexible
10	KENDA COMPLETE	For Composites, Compomers, Amalgam & Gold	Flexible
COMPOSITES			
11	KENDA MAXIMUS	For All Composites & Hybrid-Ceramics	Flexible
12	KENDA DELUXE	For All Composites, Compomers & Hybrid-Ceramics	Flexible
13	KENDA HYBRID	Optimized for Hybrid Composites	Hard
14	KENDA MICROFILL	Optimized for Microfilled Composites; Removing of Cement Remnants	Flexible
15	KENDA C.G.I.	For Universal Composites, Compomers & Glass-Ionomers	Flexible
CERAMICS			
16	KENDA ZIRCOVIS	For All Oxide & Silicate Ceramics - Optimized for Zirconia	Hard
17	KENDA ALL CERAMIC	For All Oxide & Silicate Ceramics	Flexible
18	KENDA DIAGUM	Optimized for Feldspathic Porcelain	Hard
ENAMEL			
19	KENDA PROPHY	Tooth Cleaning & Removal of Extrinsic Stains without Polishing Paste	Flexible
ACRYLICS			
20	KENDA PLANUS	For Acrylic Dentures, Retainers, Mouth Guards and Snoring Mouth Pieces	Flexible
AMALGAM AND METALS			
21	KENDA AMALGAM	For Amalgam & Precious Metals	Flexible
22	KENDA PRECIOUS METALS	For Titanium Alloys, Precious & Non-Precious Metals	Hard

CLINICAL RA - PLASTIC SHANK DISPOSABLE

PAGE	PRODUCT	APPLICATION	FLEXIBILITY
MULTIPLE MATERIALS			
23	KENDA NOBILIS	For Composites, Ceramics & Hybrid-Ceramics	Flexible
24	KENDA UNICUS	For Composites, Ceramics & Hybrid-Ceramics	Flexible
COMPOSITES			
25	KENDA HYBRID	Optimized for Hybrid Composites	Flexible
26	KENDA G.P.S.	For Composites, Compomers & Hybrid-Ceramics	Flexible
27	KENDA BIG JOB	For All Composites	Flexible
28	KENDA MICROFILL	Optimized for Microfilled Composites	Flexible
29	KENDA C.G.I.	For Universal Composites, Compomers & Glass-Ionomers	Flexible
ENAMEL			
30	KENDA PROPHY	Tooth Cleaning & Removal of Extrinsic Stains without Polishing Paste	Flexible





Grit size range specified for each step according to ISO guidelines.
Cutting properties of the same grit size vary depending on the shape and material of the abrasive particles.



The diamond symbolizes that the specific step is diamond impregnated. Diamond polishers are the most efficient, reducing polishing time and achieving the best surface results and gloss values.

PRODUCTS APPLICATION RANGE

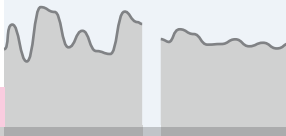CONTOURING ROUGH SURFACE		FINISHING & POLISHING SMOOTH SURFACE		ULTRA SHINE FINISH GLOSSY SURFACE		CLINICAL RA 
						PRODUCT
		32-8 μ  single step		8-4 μ  single step		MULTIPLE MATERIALS
		75-22.8 μ 1 st step		12.8-6.5 μ 2 nd step		KENDA NOBILIS
						KENDA UNICUS
						KENDA COMPLETE
						COMPOSITES
		4-1 μ  single step				KENDA MAXIMUS
		32-8 μ  1 st step		8-4 μ  2 nd step		KENDA DELUXE
		75-22.8 μ 1 st step		8-4 μ  2 nd step		KENDA HYBRID
		75-22.8 μ 1 st step		8-3 μ 2 nd step		KENDA MICROFILL
150-63 μ 1 st step		75-22.8 μ 2 nd step		8-3 μ 3 rd step		KENDA C.G.I.
						CERAMICS
		69-32 μ  1 st step		32-8 μ  2 nd step		KENDA ZIRCOVIS
		69-32 μ  1 st step		8-4 μ  2 nd step		KENDA ALL CERAMIC
		75-22.8 μ 1 st step		4-1 μ  2 nd step		KENDA DIAGUM
		75-22.8 μ single step				ENAMEL
						KENDA PROPHY
						ACRYLICS
		150-63 μ 1 st step		12.8-6.5 μ 2 nd step		KENDA PLANUS
						AMALGAM AND METALS
		75-22.8 μ 1 st step		8-3 μ 2 nd step		KENDA AMALGAM
150-63 μ 1 st step				8-4 μ  2 nd step		KENDA PRECIOUS METALS

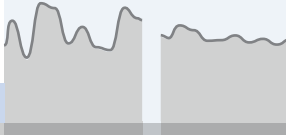CONTOURING ROUGH SURFACE		FINISHING & POLISHING SMOOTH SURFACE		ULTRA SHINE FINISH GLOSSY SURFACE		CLINICAL RA 
						PRODUCT
		32-8 μ  single step		8-4 μ  single step		MULTIPLE MATERIALS
						KENDA NOBILIS
						KENDA UNICUS
						COMPOSITES
		75-22.8 μ 1 st step		8-4 μ  2 nd step		KENDA HYBRID
		32-8 μ  1 st step		8-4 μ  2 nd step		KENDA G.P.S.
		75-22.8 μ single step				KENDA BIG JOB
		75-22.8 μ 1 st step		8-3 μ 2 nd step		KENDA MICROFILL
150-63 μ 1 st step		75-22.8 μ 2 nd step		8-3 μ 3 rd step		KENDA C.G.I.
						ENAMEL
		75-22.8 μ single step				KENDA PROPHY

Please note that the chart is only for general orientation and facilitating product choice.

Polishing results may differ depending on the specific composition of the restorative material polished.

Applied pressure and rotational speed influences the polishing result. Subject to change (product enhancements). Errors and omissions excepted.

TABLE OF CONTENTS

CLINICAL FG - STEEL SHANK MULTI-USE

PAGE	PRODUCT	APPLICATION	FLEXIBILITY
MULTIPLE MATERIALS			
31	KENDA 0306FG/0386FG	For Composites & Ceramics - On Friction Grip Shank	Hard
COMPOSITES			
31	KENDA 4006FG/5006FG	For All Composites - On Friction Grip Shank	Hard
AMALGAM AND METALS			
31	KENDA 2006FG/3006FG	For Amalgam and Precious Metals - On Friction Grip Shank	Hard

LABORATORY

PAGE	PRODUCT	APPLICATION	FLEXIBILITY
CERAMICS			
32	KENDA Z-SERIES	For All Types of Ceramics Including Zirconia, Alumina and Lithium Disilicate	Hard
33	KENDA HIGH END	For All Types of Ceramics Including Zirconia, Alumina and Lithium Disilicate	Hard
34	KENDA DIAGUM	Optimized for Feldspathic Porcelain	Hard
35	KENDA W&K 6500-7500	For Silicate Ceramics	Flexible
ACRYLICS			
36	KENDA QUEEN	For Acrylic Dentures	Flexible
37	KENDA ACRYLIC BLUE	For Acrylic Dentures	Flexible
38	KENDA RED	For Acrylic Dentures - Optimized for Partial on Chrome Cobalt Alloys	Flexible
39	KENDA YELLOW LINE	For Acrylic Dentures - Finishing Set	Flexible
40	KENDA PLANUS	For Acrylic Dentures, Retainers, Mouth Guards and Snoring Mouth Pieces	Flexible
FLEXIBLE DENTURES			
41	KENDA BESTFLEX	For All Nylon Based Flexible Dentures	Flexible
METALS			
42	KENDA PRECIOUS METALS	For All Precious, Semi and Non-Precious Metals, CoCr and Titanium	Hard
43	KENDA W&K 2500-3500	For Gold and Semi-Precious Alloys	Flexible
43	KENDA W&K 9500	For Titanium	Flexible
COMPOSITES			
44	KENDA HYBRID	For All Composite Materials	Hard
MULTIPLE MATERIALS			
45	KENDA UNIVERSAL	For Gold, Semi-Precious Metals and Composites	Hard
35	KENDA W&K 1500	For All Alloys and Acrylic Dentures	Flexible
46	KENDA 244-BLUE PASTE	For Acrylic Dentures, Flexible Partial, Composites, Gold and CoCr	

ACCESSORIES

47	KENDA PRACTICUS	Autoclavable Polisher and Bur Caddy
----	-----------------	-------------------------------------

INFORMATION

48 - 50	LABEL, INSTRUCTIONS FOR USE AND SYMBOLS
---------	---

PRODUCTS APPLICATION RANGE

CONTOURING ROUGH SURFACE



FINISHING & POLISHING SMOOTH SURFACE



ULTRA SHINE FINISH GLOSSY SURFACE



CLINICAL FG



PRODUCT

MULTIPLE MATERIALS

KENDA 0306FG/0386FG
COMPOSITES

KENDA 4006FG/5006FG
AMALGAM AND METALS

KENDA 2006FG/3006FG

CONTOURING ROUGH SURFACE



FINISHING & POLISHING SMOOTH SURFACE



ULTRA SHINE FINISH GLOSSY SURFACE



LABORATORY

PRODUCT

CERAMICS

KENDA Z-SERIES

KENDA HIGH END

KENDA DIAGUM

KENDA W&K 6500-7500

ACRYLICS

KENDA QUEEN

KENDA ACRYLIC BLUE

KENDA RED

KENDA YELLOW LINE

KENDA PLANUS

FLEXIBLE DENTURES

KENDA BESTFLEX

METALS

KENDA PRECIOUS METALS

KENDA W&K 2500-3500

KENDA W&K 9500

COMPOSITES

KENDA HYBRID

MULTIPLE MATERIALS

KENDA UNIVERSAL

KENDA W&K 1500

KENDA 244-BLUE PASTE

69-32 μ 1st step 8-4 μ 2nd step

75-22.8 μ 1st step 8-3 μ 2nd step

75-22.8 μ 1st step 8-3 μ 2nd step

164-69 μ 1st step 69-32 μ 2nd step 32-8 μ 3rd step

164-69 μ 1st step 69-32 μ 2nd step 32-8 μ 3rd step 8-4 μ 4th step

75-22.8 μ 1st step 4-1 μ 2nd step

150-63 μ 1st step 75-22.8 μ 2nd step

>150 μ 1st step 150-63 μ 2nd step

>150 μ 1st step 150-63 μ 2nd step 12.8-6.5 μ 3rd step

150-63 μ single step

150-63 μ single step

150-63 μ 1st step 12.8-6.5 μ 2nd step

150-63 μ single step

150-63 μ 1st step 8-4 μ 2nd step

75-22.8 μ 1st step 8-3 μ 2nd step

150-63 μ single step

75-22.8 μ 1st step 8-4 μ 2nd step

75-22.8 μ 1st step 8-3 μ 2nd step

150-63 μ single step

<3 μ

COMPOSITE, CERAMIC AND HYBRID-CERAMIC POLISHERS

KENDA NOBILIS



1-step diamond polishing system


One polisher only for cost-effective finishing and polishing of composite, ceramic and hybrid-ceramic restorations. The flexibility of the diamond polishers allow a gentle to teeth clean finishing and polishing. Perfect smooth surfaces will be achieved.

Recommendation:

KENDA MAXIMUS or KENDA UNICUS 1-step diamond polishers for ultimate high shine polishing.



Set REF 0300

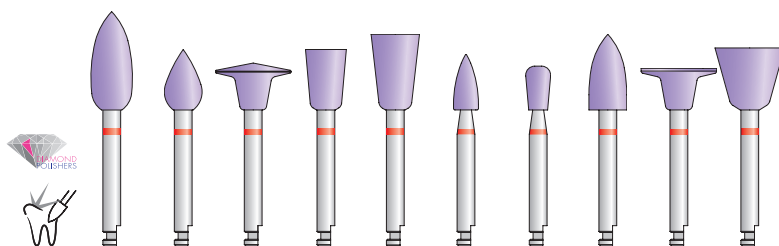

REF	reorder	0301	0302	0303	0304	0305	0306	0307	0308	0309	0310
ISO	858 204 ...	243 513 055	254 513 050	304 513 100	030 513 050	030 513 060	243 513 030	237 513 035	243 513 050	303 513 100	030 513 090
Size	Ø × L mm	5,4 × 13	5 × 8	10 × 2	5 × 8	6,5 × 10	3 × 7	3 × 6	5 × 10	10 × 1	9 × 8,5
μ	Grit	fine									
Color		violet									
rec.		7 500 - 10 000 RPM									
max.		10 000 RPM									

- Box of: 6, 12, 25, 100
- Blister of: 6 × 1, 6 × 2, 72 × 1, 72 × 2
- N Pressure: 0,3 - 2 N
- Autoclavable: 135 °C, 275 °F
- Shank: RA, stainless steel
- Medical device: class IIa

COMPOSITE, CERAMIC AND HYBRID-CERAMIC HIGH SHINE POLISHERS

KENDA UNICUS

1-step diamond polishing system


Extra-fine diamond polisher for final polishing of all composite, ceramic and hybrid-ceramic restorations. This polisher can be used on all restorations without changing the polisher. The benefits are time saving as well as cost saving.

Recommendation:

KENDA NOBILIS 1-step diamond polisher for prior finishing and polishing.



Set REF 0350

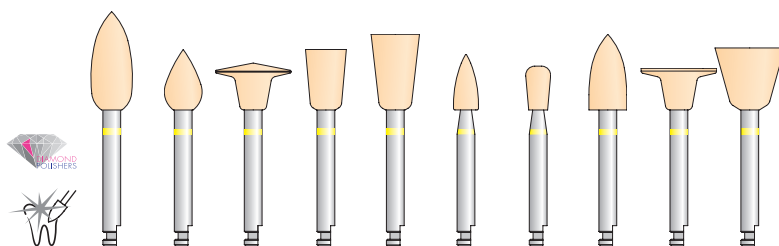

REF	reorder	0351	0352	0353	0354	0355	0356	0357	0358	0359	0360
ISO	838 204 ...	243 503 055	254 503 050	304 503 100	030 503 050	030 503 060	243 503 030	237 503 035	243 503 050	303 503 100	030 503 090
Size	Ø x L mm	5,4 x 13	5 x 8	10 x 2	5 x 8	6,5 x 10	3 x 7	3 x 6	5 x 10	10 x 1	9 x 8,5
μ	Grit	extra-fine									
	Color	orange									
	rec.	5 000 - 7 500 RPM									
	max.	10 000 RPM									

-  Box of: 6, 12, 25, 100
-  Blister of: 6 x 1, 6 x 2, 72 x 1, 72 x 2
- N** Pressure: 0,3 - 2 N
-  Autoclavable: 135 °C, 275 °F
- Shank: RA, stainless steel
- Medical device: class IIa

COMPOSITE, COMPOMER, AMALGAM & GOLD POLISHERS


KENDA COMPLETE

2-step polishing system

Universal polisher for cost-effective polishing. For all composite, compomer and amalgam restorations, for gold and dental enamel. The long-lasting prepolisher gives smooth surface finish. The high shine polisher achieves a fine surface with a long lasting high shine finish.



Set REF 0500

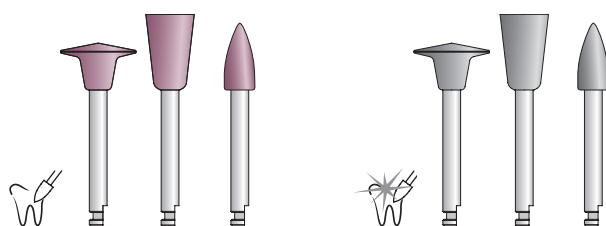

	reorder	0503	0505	0513	0553	0555	0563	
	658 204 ...	304 523 100	030 523 060	243 523 040	658 204 ...	304 513 100	030 513 060	243 513 040
	Size Ø×L mm	10×2	6,5×10	4×9	10×2	6,5×10	4×9	
	Grit	medium			fine			
	Color	light ruby colored			steel-grey			
	rec.	7 500 - 10 000 RPM			5 000 - 7 500 RPM			
	max.	10 000 RPM			10 000 RPM			

- Box of: 6, 12, 25, 100
- Blister of: 6 x 1, 6 x 2, 72 x 1, 72 x 2
- N Pressure: 0,3 - 2 N
- Autoclavable: 135 °C, 275 °F
- Shank: RA, stainless steel
- Medical device: class IIa

COMPOSITE AND HYBRID-CERAMIC POLISHERS

KENDA MAXIMUS



1-step diamond polishing system


The one step polishing technology maximises the more awkward two and three step process. Speed and contact pressure dictate the workings of KENDA MAXIMUS Polisher. By increasing the contact pressure and speed you smooth your composites and hybrid-ceramics and by decreasing the speed and pressure you create a high deep glossy shine. This one step process doesn't only give you the perfect finish but also saves you and your patient time.



Set REF 0380



«This polisher has tested to give a significant better shine and finer surface.»

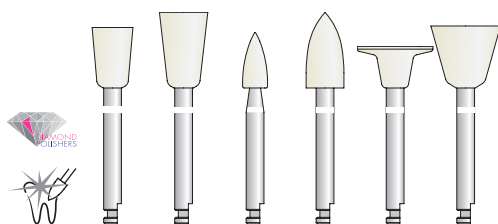

REF	reorder	0384	0385	0386	0388	0389	0390
ISO	828 204 ...	030 493 050	030 493 060	243 493 030	243 493 050	303 493 100	030 493 090
Size	Ø × L mm	5 × 8	6,5 × 10	3 × 7	5 × 10	10 × 1	9 × 8,5
μ	Grit	ultra-fine					
Color		light beige					
rec.		5 000 - 7 500 RPM					
max.		10 000 RPM					

- Box of: 6, 12, 25, 100
- Blister of: 6 × 1, 6 × 2, 72 × 1, 72 × 2
- N** Pressure: 0,3 - 2 N
- Autoclavable: 135 °C, 275 °F
- Shank: RA, stainless steel
- Medical device: class IIa

COMPOSITE, COMPOMER AND HYBRID-CERAMIC POLISHERS

KENDA DELUXE




2-step diamond polishing system

This autoclavable multiuse 2-step diamond polishing system is for extra fine smoothing and for deluxe high shine finish on all composite, compomer and hybrid-ceramic restorations.



Set REF 0400

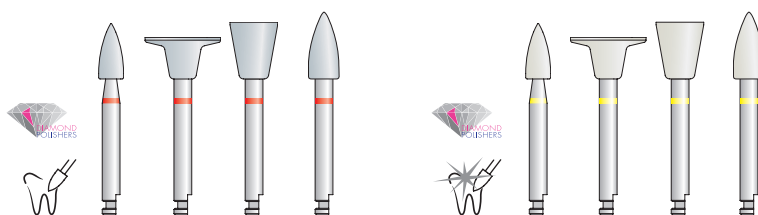

	reorder	0406	0409	0411	0413	0456	0459	0461	0463	
	858 204 ...	243 513 030	303 513 100	030 513 060	243 513 040	838 204 ...	243 503 030	303 503 100	030 503 060	243 503 040
	Size Ø x L mm	3 x 7	10 x 1	6 x 7,5	4 x 9	3 x 7	10 x 1	6 x 7,5	4 x 9	
	Grit	fine					extra-fine			
	Color	light grey					light beige			
	rec.	7 500 - 10 000 RPM					5 000 - 7 500 RPM			
	max.	10 000 RPM					10 000 RPM			

- Box of: 6, 12, 25, 100
- Blister of: 6 x 1, 6 x 2, 72 x 1, 72 x 2
- N Pressure: 0,3 - 2 N
- Autoclavable: 135 °C, 275 °F
- Shank: RA, stainless steel
- Medical device: class IIa

COMPOSITE HYBRID POLISHERS


KENDA HYBRID

2-step polishing system

Special hard polishers for all hybrid composites. The long-lasting prepolymer gives smooth surface finish. The high shine diamond polisher achieves a fine surface with a long lasting high shine finish.



Set REF 790

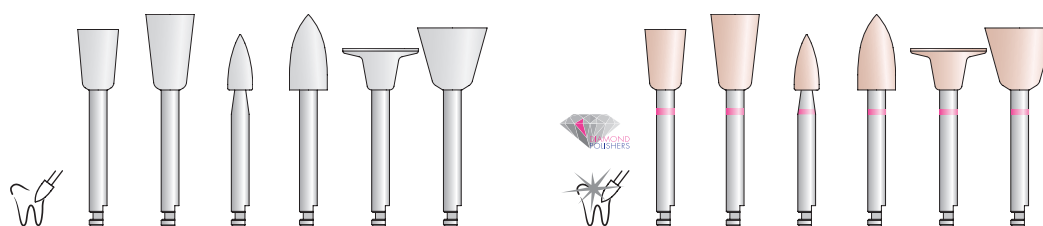

 reorder	0004	0005	0006	0008	0009	0010	0104	0105	0106	0108	0109	0110	
 632 204 ...	030 524 050	030 524 060	243 524 030	243 524 050	303 524 100	030 524 090	822 204 ...	030 494 050	030 494 060	243 494 030	243 494 050	303 494 100	030 494 090
 Size Ø×L mm	5×8	6,5×10	3×7	5×10	10×1	9×8,5	5×8	6,5×10	3×7	5×10	10×1	9×8,5	
 μ Grit	medium						ultra-fine						
 Color	whitish grey						pale pink						
 rec.	7 500 - 10 000 RPM						5 000 - 7 500 RPM						
 max.	10 000 RPM						10 000 RPM						

- Box of: 6, 12, 25, 100
- Blister of: 6×1, 6×2, 72×1, 72×2
- N Pressure: 0,3 - 2 N
- Autoclavable: 135 °C, 275 °F
- Shank: RA, stainless steel
- Medical device: class IIa

COMPOSITE MICROFILL POLISHERS

KENDA COMPOSITE MICROFILL


2-step polishing system

Top selling polisher for all microfilled composite restorations. Achieves an optimal and long lasting high shine.

Also suited for simply removing cement remnants of orthodontic bands (brackets).



Set REF 750

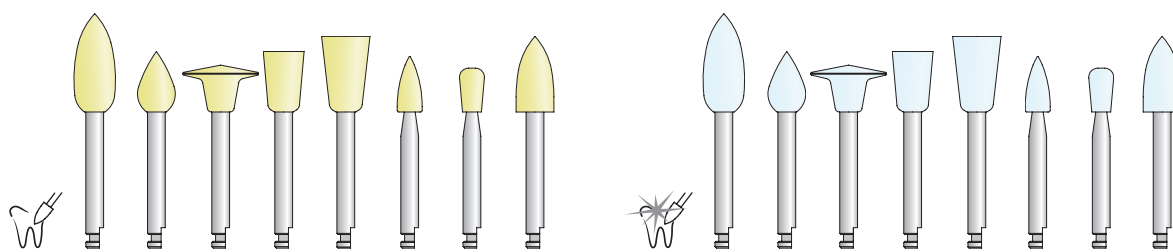

REF	reorder	4001	4002	4003	4004	4005	4006	4007	4008	5001	5002	5003	5004	5005	5006	5007	5008	
ISO	638 204 ...	243 523 055	254 523 050	304 523 100	030 523 050	030 523 060	243 523 030	237 523 035	243 523 050	638 204 ...	243 493 055	254 493 050	304 493 100	030 493 050	030 493 060	243 493 030	237 493 035	243 493 050
 Size Ø×L mm		5,4×13	5×8	10×2	5×8	6,5×10	3×7	3×6	5×10	5,4×13	5×8	10×2	5×8	6,5×10	3×7	3×6	5×10	
 Grit		medium								ultra-fine								
 Color		yellow								blue								
 rec.		7 500 - 10 000 RPM								5 000 - 7 500 RPM								
 max.		10 000 RPM								10 000 RPM								

- Box of: 6, 12, 25, 100
- Blister of: 6×1, 6×2, 72×1, 72×2
- N** Pressure: 0,3 - 2 N
- Autoclavable: 135 °C, 275 °F
- Shank: RA, stainless steel
- Medical device: class IIa

COMPOSITE & COMPOMER POLISHERS


KENDA C.G.I.


3-step polishing system

For universal composites, compomers and glass-ionomers. The C.G.I. polishing system contains specific polishers in 3 grades for contouring, finishing and polishing.



Set REF 900

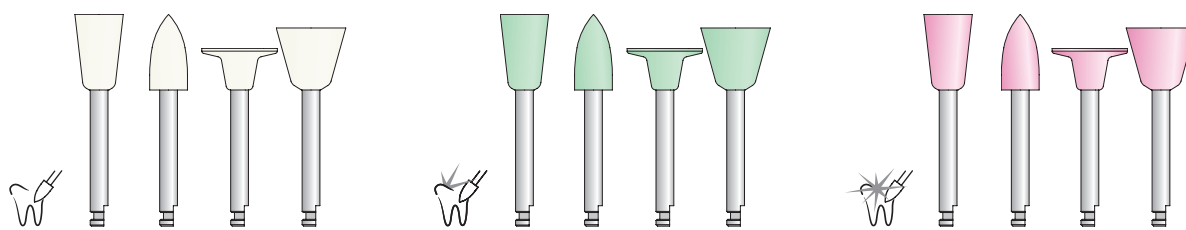

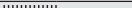REF	reorder	905-C	908-C	909-C	910-C	905-M	908-M	909-M	910-M	905-F	908-F	909-F	910-F		
ISO	638 204 ...	030 533 060	243 533 050	303 533 100	030 533 090	638 204 ...	030 523 060	243 523 050	303 523 100	030 523 090	638 204 ...	030 493 060	243 493 050	303 493 100	030 493 090
	Size Ø×L mm	6,5×10	5×10	10×1	9×8,5	6,5×10	5×10	10×1	9×8,5	6,5×10	5×10	10×1	9×8,5		
	Grit	coarse				medium				ultra-fine					
	Color	whitish white				green				pink					
	rec.	7 500 - 10 000 RPM				3 000 - 7 500 RPM				3 000 - 7 500 RPM					
	max.	10 000 RPM				10 000 RPM				10 000 RPM					

- Box of: 6, 12, 25, 100
- Blister of: 6 x 1, 6 x 2, 72 x 1, 72 x 2
- N Pressure: 0,3 - 2 N
- Autoclavable: 135 °C, 275 °F
- Shank: RA, stainless steel
- Medical device: class IIa

CERAMIC POLISHERS, OPTIMISED FOR ZIRCONIA

KENDA ZIRCOVIS




2-step diamond polishing system

Hard, robust diamond polisher, ideal for smoothing adjustments on zirconia, alumina and lithium-disilicate restorations, as well as zirconia implant components. KENDA ZIRCOVIS 2-step polishing system smooths in the first step rough surfaces, in the second step it removes fine scratches and achieves very high gloss surfaces.



Set REF 0250

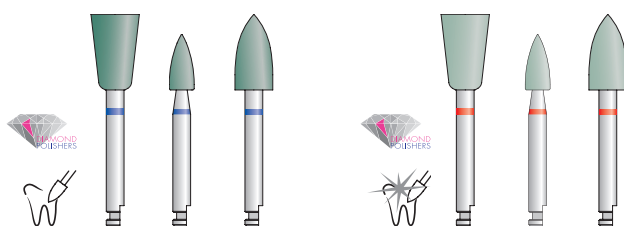

	reorder	0255	0256	0258	0265	0266	0268	
	862 204 ...	030 524 060	243 524 030	243 524 050	852 204 ...	030 514 060	243 514 030	243 514 050
	Size Ø×L mm	6,5×10	3×7	5×10	6,5×10	3×7	5×10	
	Grit	medium			fine			
	Color	dark blue-green			light blue-green			
	rec.	7 500 - 10 000 RPM			5 000 - 7 500 RPM			
	max.	10 000 RPM			10 000 RPM			

- Box of: 6, 12, 25, 100
- Blister of: 6×1, 6×2, 72×1, 72×2
- N Pressure: 0,3 - 2 N
- Autoclavable: 135 °C, 275 °F
- Shank: RA, stainless steel
- Medical device: class IIa

CERAMIC POLISHERS, FOR ALL TYPES OF CERAMICS

KENDA ALL CERAMIC




2-step diamond polishing system

Flexible polishers designed for polishing and high shine polishing of all ceramic restorations (oxide-ceramics and silicate-ceramics). Especially for polishing the margins and chairside occlusal corrections.



Set REF 0200

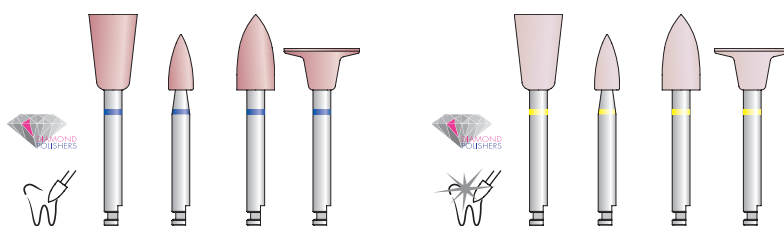

REF	reorder	0215	0216	0218	0219	0225	0226	0228	0229	
ISO	868 204 ...	030 523 060	243 523 030	243 523 050	303 523 100	838 204 ...	030 503 060	243 503 030	243 503 050	303 503 100
	2 3	Size Ø×L mm	6,5×10	3×7	5×10	10×1	6,5×10	3×7	5×10	10×1
μ	Grit	medium				extra-fine				
ℓ	Color	medium bordeau				light bordeau				
↻	rec.	7 500 - 10 000 RPM				5 000 - 7 500 RPM				
↻	max.	10 000 RPM				10 000 RPM				

- Box of: 6, 12, 25, 100
- Blister of: 6 x 1, 6 x 2, 72 x 1, 72 x 2
- N Pressure: 0,3 - 2 N
- Autoclavable: 135 °C, 275 °F
- Shank: RA, stainless steel
- Medical device: class IIa

CERAMIC POLISHERS, OPTIMISED FOR FELDSPATHIC PORCELAIN

KENDA DIAGUM


2-step polishing system

Hard polishers for excellent results on feldspathic porcelain in all forms.

Silicon carbide prepolishers for smoothing and fine diamond high shine polishers for brilliant lustre.



Set REF 780

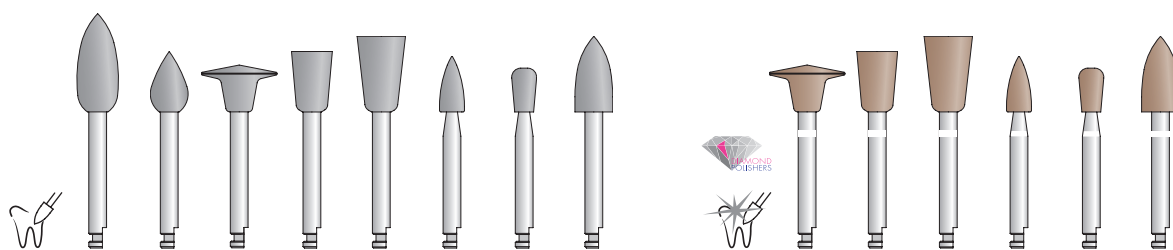

REF	reorder	6001	6002	6003	6004	6005	6006	6007	6008	X003	X004	X005	X006	X007	X008	
ISO	652 204 ...	243 534 055	254 534 050	304 534 100	030 534 050	030 534 060	243 534 030	237 534 035	243 534 050	822 204 ...	304 494 100	030 494 050	030 494 060	243 494 030	237 494 035	243 494 050
	Size Ø×L mm	5,4×13	5×8	10×2	5×8	6,5×10	3×7	3×6	5×10	10×2	5×8	6,5×10	3×7	3×6	5×10	
	Grit	coarse									ultra-fine					
	Color	grey									dark rosa					
	rec.	7 500 - 10 000 RPM									5 000 - 7 500 RPM					
	max.	10 000 RPM									10 000 RPM					

- Box of: 6, 12, 25, 100
- Blister of: 6 x 1, 6 x 2, 72 x 1, 72 x 2
- N Pressure: 0,3 - 2 N
- Autoclavable: 135 °C, 275 °F
- Shank: RA, stainless steel
- Medical device: class IIa

PROPHYLAXIS POLISHERS


KENDA PROPHY

1-step polishing system

For tooth cleaning without the use of polishing paste. Medium grit polishers for removal of extrinsic stains, caused by chromogenic agents or surface staining (plaque and calculus). To be used with gentle pressure.



Set REF 5100

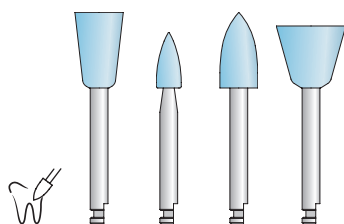

REF reorder	5105	5106	5108	5110
ISO 658 204 ...	030 523 060	243 523 030	243 523 050	030 523 090
Size Ø×L mm	6,5×10	3×7	5×10	9×8,5
μ Grit	medium			
Color	blue			
rec.	7 500 - 10 000 RPM			
max.	10 000 RPM			

- Box of: 6, 12, 25, 100
- Blister of: 6×1, 6×2, 72×1, 72×2
- N** Pressure: 0,3 - 2 N
- Autoclavable: 135 °C, 275 °F
- Shank: RA, stainless steel
- Medical device: class IIa

ACRYLIC DENTURE POLISHERS

KENDA PLANUS


2-step smoothing and polishing system

This 2-step autoclavable polisher is designed for chairside smoothing of pressure spots on acrylic dentures and retainers. Also for work on mouth guards (bruxism treatment) and snoring mouth pieces.

Available in RA and HP version (see page 40).



Set REF 0600

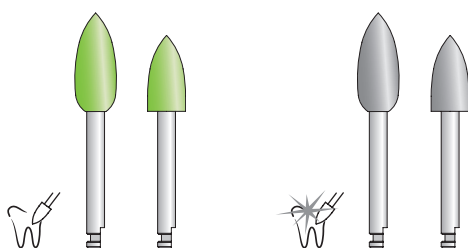

REF	reorder	0601	0608	0611	0618
ISO	638 204 ...	243 533 055	243 533 050	638 204 ...	243 513 055
Size Ø × L mm		5,4 × 13	5 × 10	5,4 × 13	5 × 10
μ Grit		coarse		fine	
Color		lime-green		grey	
rec.		7 500 - 10 000 RPM		5 000 - 7 500 RPM	
max.		10 000 RPM		10 000 RPM	

- Box of: 6, 12, 25, 100
- Pressure: 0,3 - 2 N
- Autoclavable: 135 °C, 275 °F
- Shank: RA, stainless steel
- Medical device: class IIa

AMALGAM POLISHERS

KENDA AMALGAM


2-step polishing system

Long lasting silicone polisher with fine adjustments of integrated polish care for fast and brilliant polishing of amalgam and precious metals.

The prepolisher achieves a smooth surface and the high shine polisher a perfect mirror lustre.



Set REF 700

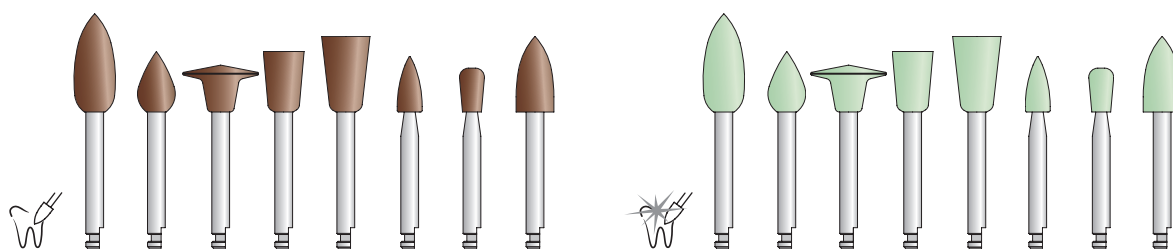

REF	reorder	2001	2002	2003	2004	2005	2006	2007	2008	3001	3002	3003	3004	3005	3006	3007	3008	
ISO	658 204 ...	243 523 055	254 523 050	304 523 100	030 523 050	030 523 060	243 523 030	237 523 035	243 523 050	658 204 ...	243 513 055	254 513 050	304 513 100	030 513 050	030 513 060	243 513 030	237 513 035	243 513 050
 Size Ø×L mm		5,4×13	5×8	10×2	5×8	6,5×10	3×7	3×6	5×10	5,4×13	5×8	10×2	5×8	6,5×10	3×7	3×6	5×10	
 Grit		medium								fine								
 Color		brown								green								
 rec.		7 500 - 10 000 RPM								5 000 - 7 500 RPM								
 max.		10 000 RPM								10 000 RPM								

- Box of: 6, 12, 25, 100
- Blister of: 6×1, 6×2, 72×1, 72×2
- N** Pressure: 0,3 - 2 N
- Autoclavable: 135 °C, 275 °F
- Shank: RA, stainless steel
- Medical device: class IIa

METAL POLISHERS

KENDA PRECIOUS METALS

2-step polishing system

For finishing, pre-polishing and high shine. Results in a matter of seconds.

The high shine polisher has microfine diamond particles with self sharpening cutting edges for the highest efficiency on precious and non-precious metals as well as titanium alloys.



Set REF 740



	reorder	1004	1005	1006	1008	1009	1010	1104	1105	1106	1108	1109	1110	
	652 204 ...	030 524 050	030 524 060	243 524 030	243 524 050	303 524 100	030 524 090	832 204 ...	030 504 050	030 504 060	243 504 030	243 504 050	303 504 100	030 504 090
	Size Ø×L mm	5×8	6,5×10	3×7	5×10	10×1	9×8,5	5×8	6,5×10	3×7	5×10	10×1	9×8,5	
	Grit	medium						extra-fine						
	Color	black						green						
	rec.	7 500 - 10 000 RPM						5 000 - 7 500 RPM						
	max.	10 000 RPM						10 000 RPM						

- Box of: 6, 12, 25, 100
- Blister of: 6×1, 6×2, 72×1, 72×2
- N Pressure: 0,3 - 2 N
- Autoclavable: 135 °C, 275 °F
- Shank: RA, stainless steel
- Medical device: class IIa

COMPOSITE, CERAMIC AND HYBRID-CERAMIC POLISHERS

DISPOSABLE KENDA NOBILIS

1-step diamond polishing system


For single patient use.

One polisher only for cost-effective finishing and polishing of composite, ceramic and hybrid-ceramic restorations. The flexibility of the diamond polishers allow a gentle to teeth clean finishing and polishing. Perfect smooth surfaces will be achieved.

Disposable polisher saves sterilisation cost and time.



Set REF PS0300

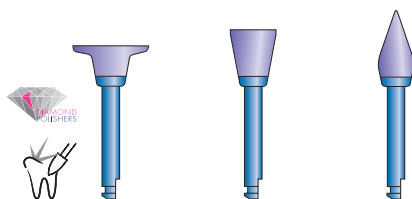

REF	reorder	PS0309	PS0311	PS0312
ISO	858 242 ...	303 513 100	030 513 060	243 513 045
Size Ø x L mm		10 x 1	6 x 7,5	5 x 10
μ Grit		fine		
Color		violet		
rec.		5 000 RPM		
max.		5 000 RPM		

-  Box of: 6, 12, 25, 100
- N** Pressure: 0,3 - 1 N
- Shank: RA, medical plastic
- Medical device: class IIa
-  disposable

COMPOSITE, CERAMIC AND HYBRID-CERAMIC HIGH SHINE POLISHERS

DISPOSABLE KENDA UNICUS

1-step diamond polishing system

For single patient use.


Extra-fine diamond polisher for final polishing of all composite, ceramic and hybrid-ceramic restorations. This polisher can be used on all restorations without changing the polisher. The benefits are time saving as well as cost saving. Disposable polisher saves sterilisation cost and time.

Recommendation:

Disposable KENDA NOBILIS 1-step diamond polisher for finishing and polishing.



Set REF PS0350

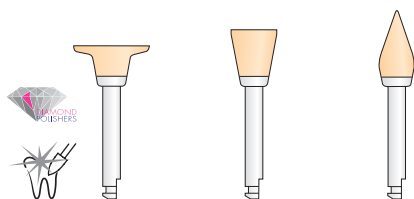

REF	reorder	PS0359	PS0361	PS0362
ISO	838 242 ...	303 503 100	030 503 060	243 503 045
Size Ø x L mm		10 x 1	6 x 7,5	5 x 10
μ Grit		extra-fine		
Color		orange		
rec.		5 000 RPM		
max.		5 000 RPM		

-  Box of: 6, 12, 25, 100
- N** Pressure: 0,3 - 1 N
- Shank: RA, medical plastic
- Medical device: class IIa
-  disposable

COMPOSITE HYBRID POLISHERS

DISPOSABLE KENDA HYBRID


2-step polishing system

For single patient use.

Special polishers for all hybrid composites. The long-lasting prepolymer gives smooth surface finish. The high shine diamond polisher achieves a fine surface with a long lasting high shine finish. Disposable polisher saves sterilisation cost and time.



Set REF PS790

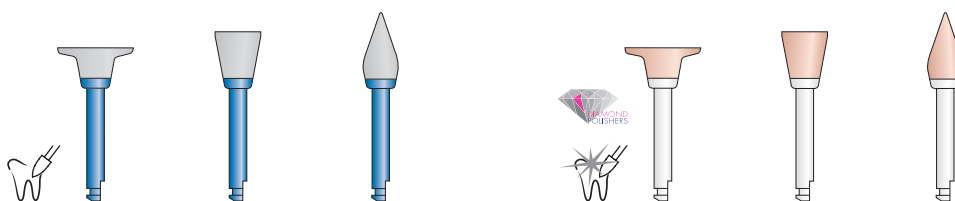

	reorder	PS0009	PS0011	PS0012	PS0109	PS0111	PS0112	
	638 242 ...	303 523 100	030 523 060	243 523 045	828 242 ...	303 493 100	030 493 060	243 493 045
	Size Ø×L mm	10×1	6×7,5	5×10	10×1	6×7,5	5×10	
	Grit	medium			ultra-fine			
	Color	whitish grey			pale pink			
	rec.	5 000 RPM			5 000 RPM			
	max.	5 000 RPM			5 000 RPM			

- Box of: 6, 12, 25, 100
- Pressure: 0,3 - 1 N
- Shank: RA, medical plastic
- Medical device: class IIa
- disposable

COMPOSITE, COMPOMER AND HYBRID-CERAMIC POLISHERS

DISPOSABLE KENDA G.P.S.

2-step diamond polishing system

For single patient use.


For all composite, compomer and hybrid-ceramic restorations. Prepolariser for smoothing and final polisher for finest surfaces and highest long lasting shine.

Disposable polisher saves sterilisation cost and time.



Set REF PS0400

**THE DENTAL
ADVISOR**
++++ ½

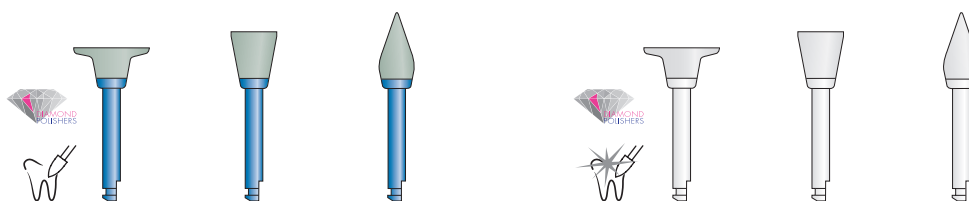

	reorder	PS0409	PS0411	PS0412	PS0459	PS0461	PS0462	
	858 242 ...	303 513 100	030 513 060	243 513 045	838 242 ...	303 503 100	030 503 060	243 503 045
	Size Ø×L mm	10×1	6×7,5	5×10	10×1	6×7,5	5×10	
	Grit	fine			extra-fine			
	Color	light green			light beige			
	rec.	5 000 RPM			5 000 RPM			
	max.	5 000 RPM			5 000 RPM			

-  Box of: 6, 12, 25, 100
- N** Pressure: 0,3 - 1 N
- Shank: RA, medical plastic
- Medical device: class IIa
-  disposable

Edition 08/2016

FINISHING AND POLISHING OF ALL COMPOSITES

DISPOSABLE KENDA BIG JOB

1-step polishing system

For single patient use.


This fast working polisher is designed for large job multi restorations. Ideal for fine surface polishing of all composite restorations. Disposable polisher saves sterilisation cost and time.

Recommendation:

For ultimate shine use G.P.S. diamond hi-shine polisher.



Set REF PS4100

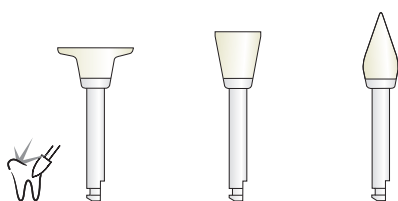

REF	reorder	PS4109	PS4111	PS4112
ISO	638 242 ...	303 523 100	030 523 060	243 523 045
Size Ø × L mm		10 × 1	6 × 7,5	5 × 10
μ Grit		medium		
Color		white		
rec.		5 000 RPM		
max.		5 000 RPM		

- Box of: 6, 12, 25, 100
- N** Pressure: 0,3 - 1 N
- Shank: RA, medical plastic
- Medical device: class IIa
- disposable

COMPOSITE MICROFILL POLISHERS

DISPOSABLE KENDA COMPOSITE MICROFILL

2-step polishing system


For single patient use.

Top selling polisher for all microfilled composite restorations. Achieves an optimal and long lasting high shine. Also suited for simply removing cement remnants of orthodontic bands (brackets).

Disposable polisher saves sterilisation cost and time.



Set REF PS750

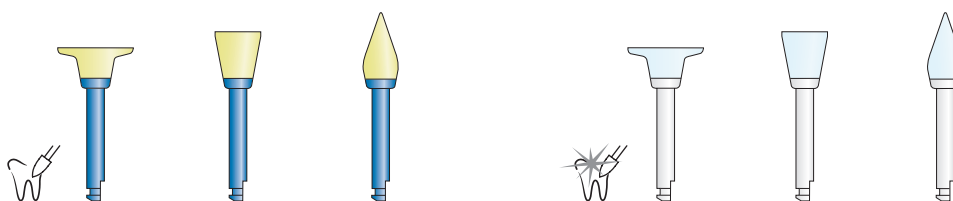

	reorder	PS4009	PS4011	PS4012	PS5009	PS5011	PS5012	
	638 242 ...	303 523 100	030 523 060	243 523 045	638 242 ...	303 493 100	030 493 060	243 493 045
	Size Ø×L mm	10×1	6×7,5	5×10	10×1	6×7,5	5×10	
	Grit	medium			ultra-fine			
	Color	yellow			blue			
	rec.	5 000 RPM			5 000 RPM			
	max.	5 000 RPM			5 000 RPM			

- Box of: 6, 12, 25, 100
- N** Pressure: 0,3 - 1 N
- Shank: RA, medical plastic
- Medical device: class IIa
- disposable

COMPOSITE & COMPOMER POLISHERS

DISPOSABLE KENDA C.G.I.

3-step polishing system




For single patient use.

For universal composites, compomers and glassionomers. The C.G.I. polishing system contains specific polishers in 3 grades for contouring, finishing and polishing.

Disposable polisher saves sterilisation cost and time.



Set REF PS900

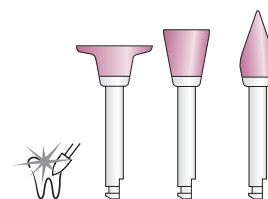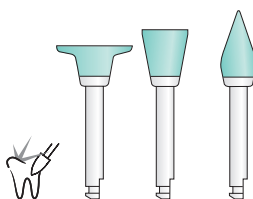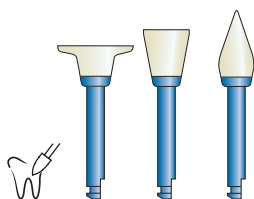

REF	reorder	PS909-C	PS911-C	PS912-C	PS909-M	PS911-M	PS912-M	PS909-F	PS911-F	PS912-F		
ISO	638 242 ...	303 533 100	030 533 060	243 533 045	638 242 ...	303 523 100	030 523 060	243 523 045	638 242 ...	303 493 100	030 493 060	243 493 045
 Size Ø×L mm		10×1	6×7,5	5×10	10×1	6×7,5	5×10	10×1	6×7,5	5×10		
 Grit		coarse			medium			ultra-fine				
 Color		whitish white			green			pink				
 rec.		5 000 RPM			5 000 RPM			3 000 RPM				
 max.		5 000 RPM			5 000 RPM			5 000 RPM				

- Box of: 6, 12, 25, 100
- N** Pressure: 0,3 - 1 N
- Shank: RA, medical plastic
- Medical device: class IIa
- disposable

PROPHYLAXIS POLISHERS

DISPOSABLE KENDA PROPHY

1-step polishing system


For single patient use.

For tooth cleaning without the use of polishing paste. Medium grit polishers for removal of extrinsic stains, caused by chromogenic agents or surface staining (plaque and calculus). To be used with gentle pressure.

Disposable polisher saves sterilisation cost and time.



Set REF PS5100

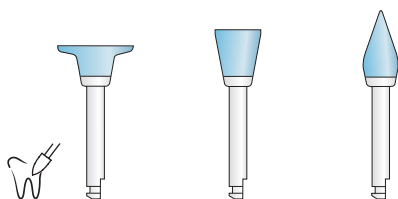


REF	reorder	PS5109	PS5111	PS5112
ISO	658 242 ...	303 523 100	030 523 060	243 523 045
Size	Ø × L mm	10 × 1	6 × 7,5	5 × 10
μ	Grit	medium		
Color		blue		
rec.		5 000 RPM		
max.		5 000 RPM		

- Box of: 6, 12, 25, 100
- N** Pressure: 0,3 - 1 N
- Shank: RA, medical plastic
- Medical device: class IIa
- Ⓢ disposable


FRICION GRIP POLISHING SYSTEMS

KENDA FG SERIES

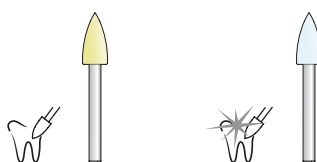


2-step polishing system for all composites

REF	reorder	4006FG	5006FG
ISO	632 313 ...	243 524 030	632 313 ... 243 494 030
Size Ø × L mm	3 × 7	3 × 7	3 × 7
μ Grit	medium	ultra-fine	
Color	yellow	blue	
rec.	20 000 RPM	10 000 RPM	
max.	20 000 RPM	20 000 RPM	

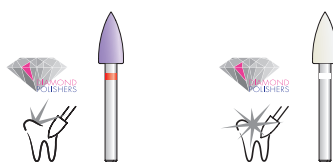

2-step diamond polishing system for all composites and ceramics. For very efficient polishing.

REF	reorder	0306FG	0386FG
ISO	852 313 ...	243 514 030	822 313 ... 243 494 030
Size Ø × L mm	3 × 7	3 × 7	3 × 7
μ Grit	fine	ultra-fine	
Color	violet	light beige	
rec.	20 000 RPM	10 000 RPM	
max.	20 000 RPM	20 000 RPM	

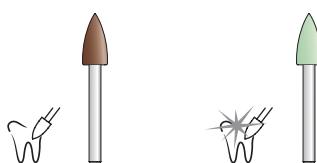

2-step polishing system for amalgam and precious metals.

REF	reorder	2006FG	3006FG
ISO	652 313 ...	243 524 030	652 313 ... 243 514 030
Size Ø × L mm	3 × 7	3 × 7	3 × 7
μ Grit	medium	fine	
Color	brown	green	
rec.	20 000 RPM	10 000 RPM	
max.	20 000 RPM	20 000 RPM	

- Box of: 5
- Pressure: 0,3 - 2 N
- Autoclavable: 135 °C, 275 °F
- Shank: FG, stainless steel
- Medical device: class IIa

CERAMIC POLISHERS, FOR ALL TYPES OF CERAMICS

KENDA Z-SERIES




3-step diamond polishing system

For best results on high performance ceramic, zirconia, alumina and lithium disilicate.

First step allows fast contouring and trimming, followed by the second step smoothing and polishing and finally with the third step for high shine finishing.



Set REF Z-000

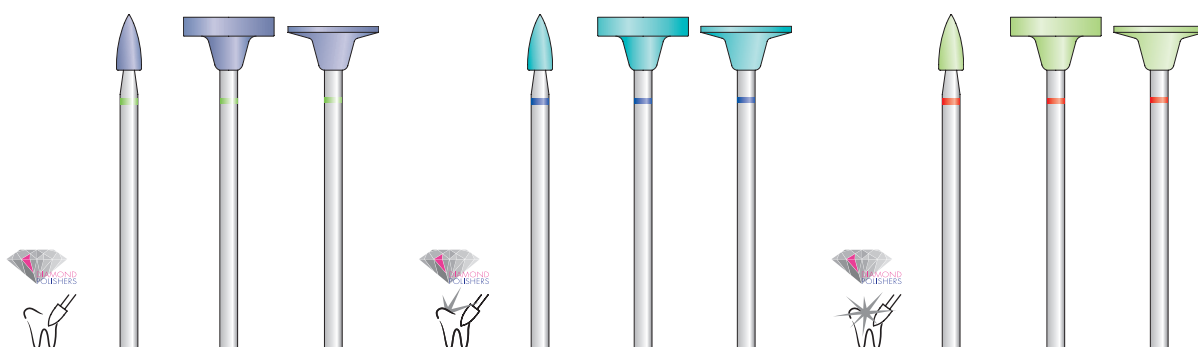

REF	reorder	Z-006-C	Z-010-C	Z-014-C	Z-006-M	Z-010-M	Z-014-M	Z-006-F	Z-010-F	Z-014-F
ISO	872 104 ...	243 534 030	041 534 120	041 534 120	862 104 ...	243 524 030	041 524 100	852 104 ...	243 514 030	041 514 120
Size	Ø x L mm	3 x 7	12 x 2,5	12 x 2	3 x 7	12 x 2,5	12 x 2	3 x 7	12 x 2,5	12 x 2
μ	Grit	coarse			medium			fine		
	Color	violet			petrol			lime		
rec.		7 500 - 10 000 RPM			7 500 - 10 000 RPM			5 000 - 7 500 RPM		
max.		10 000 RPM			10 000 RPM			10 000 RPM		

Box of: 3
Pressure: 0,3 - 2 N
Shank: HP, stainless steel

CERAMIC POLISHERS, FOR ALL TYPES OF CERAMICS

KENDA HIGH END


4-step diamond polishing system

The power disc represents a new dimension for the complete treatment of all types of ceramics, especially pressed ceramic and zirconia trimming, smoothing, finishing and high shine polishing.

The special geometry of the disc allows better access to the interdental areas. The extra coarse disc is ideal for fast trimming and smoothing of pressed ceramic and zirconia.



Set REF 626

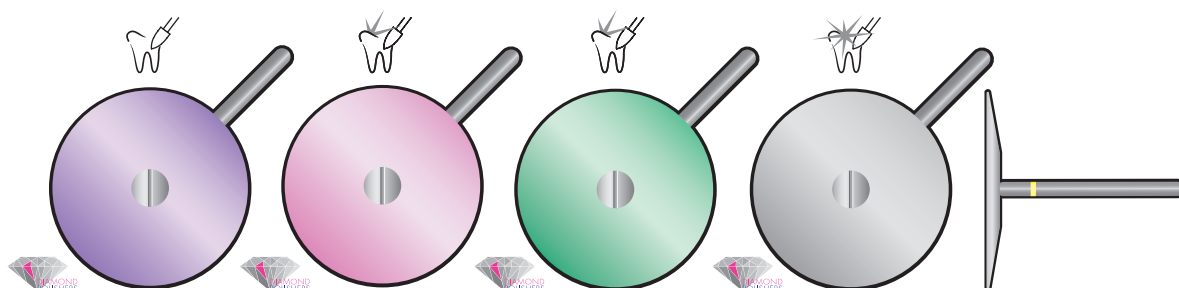

REF	reorder	626-C	626-M	626-F	626-XF
ISO		872 104 303 534 260	862 104 303 524 260	852 104 303 514 260	832 104 303 504 260
Size Ø × L mm		26 × 2	26 × 2	26 × 2	26 × 2
μ Grit		coarse	medium	fine	extra-fine
Color		violet	pink	green	grey
rec.		7 500 - 10 000 RPM	7 500 - 10 000 RPM	7 500 - 10 000 RPM	5 000 - 7 500 RPM
max.		10 000 RPM	10 000 RPM	10 000 RPM	10 000 RPM

Box of: 1
N Pressure: 0,3 - 2 N
 Shank: HP, stainless steel

CERAMIC POLISHERS, OPTIMISED FOR FELDSPATHIC PORCELAIN

KENDA DIAGUM, HP


2-step polishing system

2-step polishing system for ceramics.

Grey pre-polisher for contouring and smoothing. Dark rosa hi-shine polisher to achieve rapidly a high lustre and natural aesthetics due to its content of ultra-fine diamonds.



Set REF 781

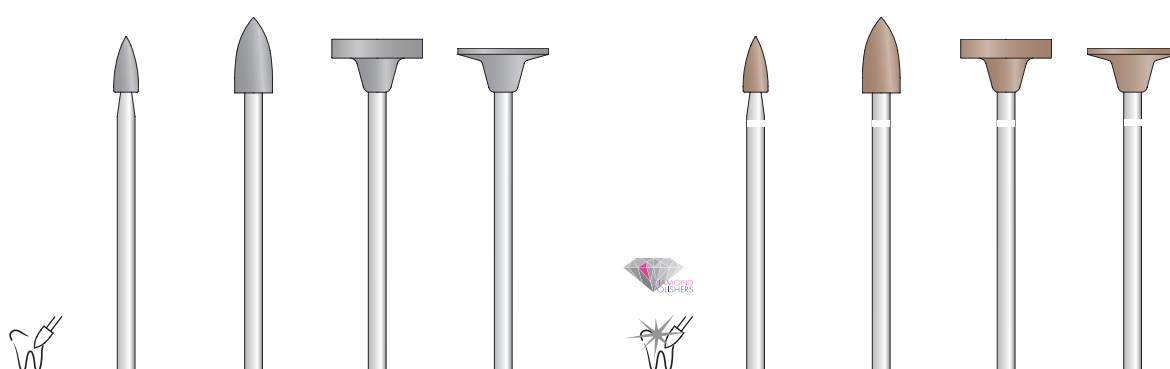


REF	reorder	6106HP	6108HP	6110HP	6114HP		X006HP	X008HP	X010HP	X014HP
ISO	652 104 ...	243 534 030	243 534 050	041 534 120	041 534 120	822 104 ...	243 494 030	243 494 050	041 494 120	041 494 120
Size	Ø x L mm	3 x 7	5 x 10	12 x 2,5	12 x 2		3 x 7	5 x 10	12 x 2,5	12 x 2
μ	Grit	coarse					ultra fine			
Color		grey					dark rosa			
rec.		7 500 - 10 000 RPM					5 000 - 7 500 RPM			
max.		10 000 RPM					10 000 RPM			

Box of: 6, 12, 25, 100
N Pressure: 0,3 - 2 N
 Shank: HP, stainless steel

CERAMIC / UNIVERSAL POLISHERS UNMOUNTED


KENDA WHEEL & KNIFE

Ceramic prepolisher for fast smoothing of all silicate ceramics.

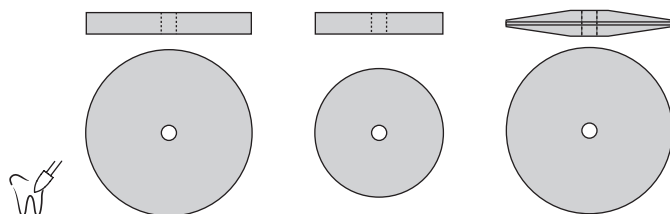


REF reorder	6522-R	6517-R	6522-L
ISO 658 900 ...	372 533 220	372 533 170	303 533 220
Size Ø×L mm	22×3	17×3	22×3
Grit	coarse		
Color	light grey		
rec.	20 000 RPM		
max.	20 000 RPM		

Ceramic high shine polisher for final shine.

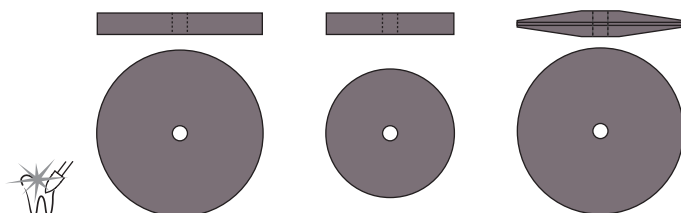

REF reorder	7522-R	7517-R	7522-L
ISO 658 900 ...	372 513 220	372 513 170	303 513 220
Size Ø×L mm	22×3	17×3	22×3
Grit	fine		
Color	pale pink		
rec.	10 000 RPM		
max.	20 000 RPM		

Universal polisher, for polishing all alloys and acrylic dentures.

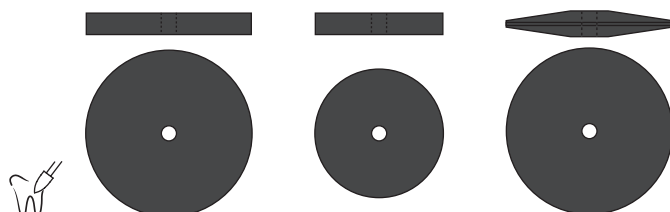

REF reorder	1522-R	1517-R	1522-L
ISO 658 900 ...	372 523 220	372 523 170	303 523 220
Size Ø×L mm	22×3	17×3	22×3
Grit	medium		
Color	black		
rec.	20 000 RPM		
max.	20 000 RPM		

Box of: 10, 100
N Pressure: 0,3 - 2 N


ACRYLIC DENTURE POLISHERS

KENDA QUEEN


2-step finishing and polishing system

Green: Ideal for contouring of partial and full dentures as well as orthodontic appliances.

Black: Ideal for fast smoothing and polishing of partial and full dentures as well as orthodontic appliances.

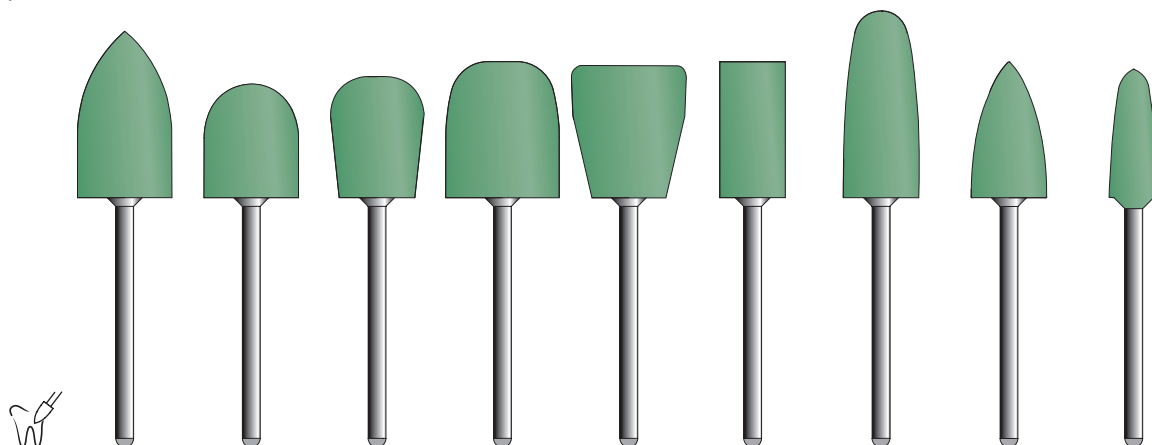

REF	reorder	9008	9009	9010	9011	9012	9013	9014	9015	9016
ISO	658 104 ...	245 533 125	144 533 125	237 533 120	201 533 150	012 533 150	114 533 085	204 533 100	243 533 100	213 533 065
Size	Ø × L mm	12,5 × 22	12,5 × 15	12 × 17	15 × 18	15 × 18	8,5 × 18	10 × 24,5	10 × 18	6,5 × 17
μ	Grit	coarse								
	Color	green								
rec.		10 000 - 15 000 RPM								
max.		15 000 RPM								

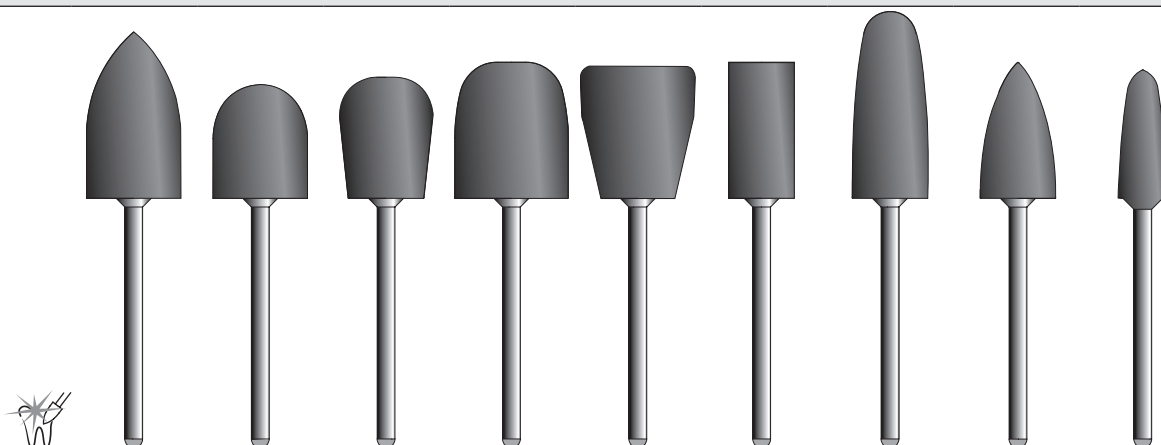

REF	reorder	7008	7009	7010	7011	7012	7013	7014	7015	7016
ISO	658 104 ...	245 523 125	144 523 125	237 523 120	201 523 150	012 523 150	114 523 085	204 523 100	243 523 100	213 523 065
Size	Ø × L mm	12,5 × 22	12,5 × 15	12 × 17	15 × 18	15 × 18	8,5 × 18	10 × 24,5	10 × 18	6,5 × 17
μ	Grit	medium								
	Color	black								
rec.		10 000 - 15 000 RPM								
max.		15 000 RPM								


Box of: 6, 12, 25, 100
N Pressure: 0,3 - 2 N
 Shank: HP, stainless steel

ACRYLIC DENTURE POLISHERS


KENDA ACRYLIC BLUE

3-step polishing system

Blue: ideal for fast trimming of acrylics
Light blue: ideal for fast smoothing of acrylics
Ochre: ideal for a high gloss finish

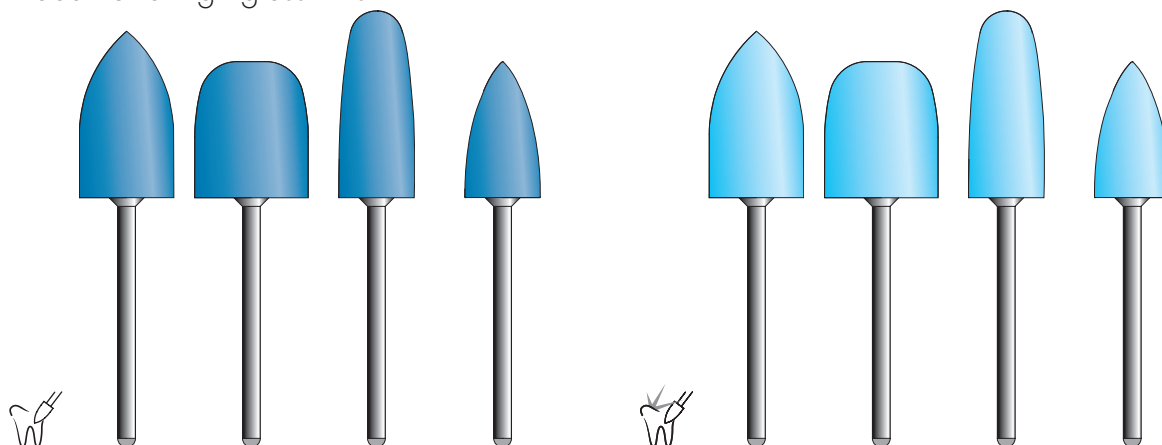

 reorder	8308	8311	8314	8315	8408	8411	8414	8415	
 638 104 ...	245 543 125	201 543 150	204 543 100	243 543 100	638 104 ...	245 523 125	201 523 150	204 523 100	243 523 100
 Size Ø×L mm	12,5×22	15×18	10×24,5	10×18		12,5×22	15×18	10×24,5	10×18
 Grit	very coarse				medium				
 Color	blue				light blue				
 rec.	10 000 - 15 000 RPM				10 000 - 15 000 RPM				
 max.	15 000 RPM				15 000 RPM				

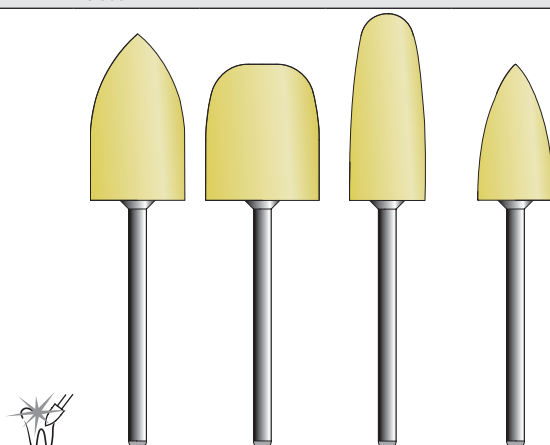

REF	reorder	8508	8511	8514	8515
ISO	638 104 ...	245 503 125	201 503 150	204 503 100	243 503 100
Size	Ø × L mm	12,5 × 22	15 × 18	10 × 24,5	10 × 18
μ	Grit	extra-fine			
	Color	ochre			
rec.		5 000 - 10 000 RPM			
max.		15 000 RPM			

Box of: 6, 12, 25, 100
Pressure: 0,3 - 2 N
Shank: HP, stainless steel

ACRYLIC DENTURE POLISHERS

KENDA RED

1-step polishing system


For adjusting and polishing dentures and flexible partials on chrome cobalt alloys without damaging the alloy, mounted and unmounted.

Recommendation:

For high shine use KENDA 244-BLUE polishing paste and SLIMLINE brush.



REF RED14.006

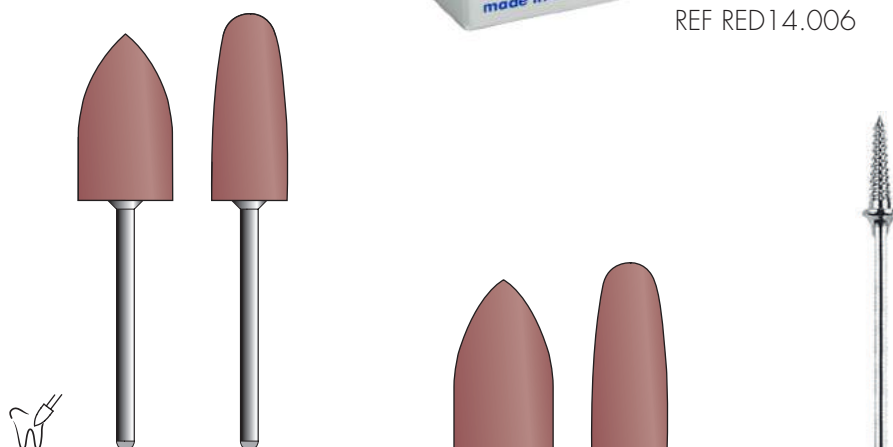

REF	reorder	RED 8	RED 14	10008	10014	301-L
ISO	658 104 ...	245 523 125	204 523 100	658 900 ...	245 523 125	204 523 100
Size	Ø × L mm	12,5 × 22	10 × 24,5	12,5 × 22	10 × 24,5	5 × 13
μ	Grit	medium		medium		
	Color	red		red		
rec.		10 000 - 15 000 RPM		10 000 RPM		
max.		15 000 RPM		15 000 RPM		

Box of: 6, 12, 25, 100

N Pressure: 0,3 - 2 N

Shank: HP, stainless steel

ACRYLIC DENTURE FINISHING SET

KENDA YELLOW LINE


Y-014: ideal for eliminating pressure spots fast and efficient.

Y-018: the rounded edges of the instrument make it possible to smooth without grooves, welts and crusts.

Y-019: super flexible brush – disc for the cleaning and polishing of the interdental areas.



Set REF Y-000

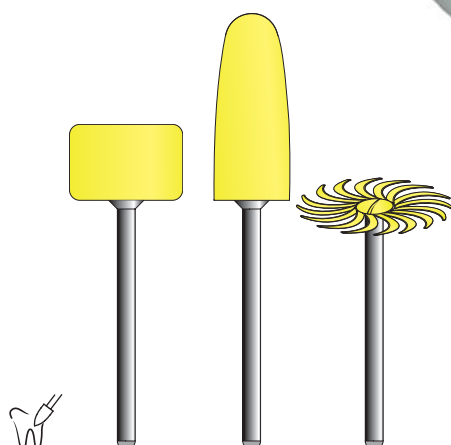

REF	reorder	Y-018	Y-014	Y-019
ISO	638 104 ...	056 523 150	204 523 100	543 523 190
Size Ø × L mm		15 × 10	10 × 24,5	19 × 0,5
μ Grit		medium		
Color		yellow		
rec.		7 500 - 10 000 RPM		
max.		10 000 RPM		

Box of: 3
N Pressure: 0,3 - 2 N
 Shank: HP, stainless steel

ACRYLIC DENTURE POLISHERS

KENDA PLANUS, HP


2-step smoothing and polishing system

This 2-step autoclavable polisher are designed for smoothing of pressure spots on acrylic dentures and retainers. Also for work on mouth guards (bruxism treatment) and snoring mouth pieces. The special medical silicone of this HP polisher allows autoclaving at 135 °C when used chairside.

Available in RA (see page 20) and HP version.



Set REF 6700

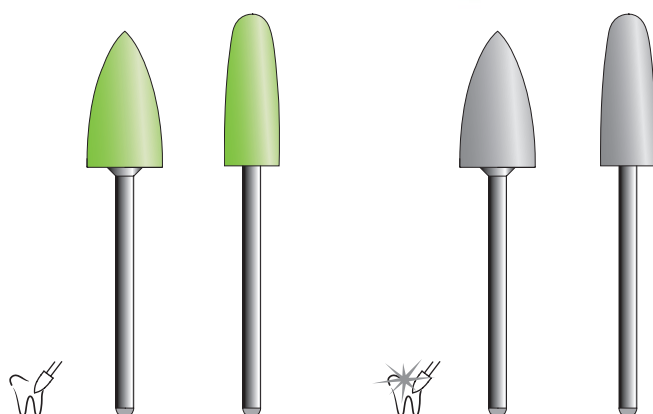

REF	reorder	6715	6720	6715	6720
ISO	638 104 ...	243 533 100	497 533 070	638 104 ...	243 513 100
Size Ø x L mm		10 x 18	10 x 18	10 x 18	10 x 18
Grit		coarse		fine	
Color		lime-green		grey	
rec.		10 000 - 15 000 RPM		10 000 - 15 000 RPM	
max.		15 000 RPM		15 000 RPM	
Box of:	6, 12, 25, 100	Autoclavable:		135 °C, 275 °F	
Pressure:	0,3 - 2 N	Shank:		HP, stainless steel	

FLEXIBLE DENTURE POLISHERS

KENDA BESTFLEX

1-step polishing system


For polishing of all nylon based flexible dentures.

Recommendation:

For high shine use KENDA 244-BLUE polishing paste and SLIMLINE brush.



Set REF 6500

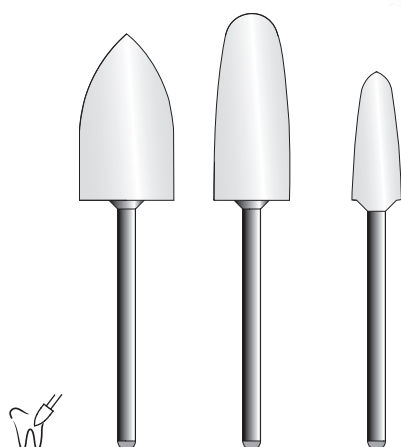

REF	reorder	6508	6514	6516
ISO	658 104 ...	245 533 125	204 533 100	213 533 065
Size	Ø × L mm	12,5 × 22	10 × 24,5	6,5 × 17
μ	Grit	coarse		
Color		white grey		
rec.		10 000 - 15 000 RPM		
max.		15 000 RPM		
Box of:		6, 12, 25, 100		
N Pressure:		0,3 - 2 N		
Shank:		HP, stainless steel		

METAL POLISHERS


KENDA PRECIOUS UNIVERSAL METALS, HP

2-step polishing system

2-step polishing system for all alloys of precious, semiprecious and non precious metals, CoCr and titanium alloys, HP. Black pre-prolisher for smoothing and the green diamond polisher for highest gloss in a matter of seconds.



Set REF 741

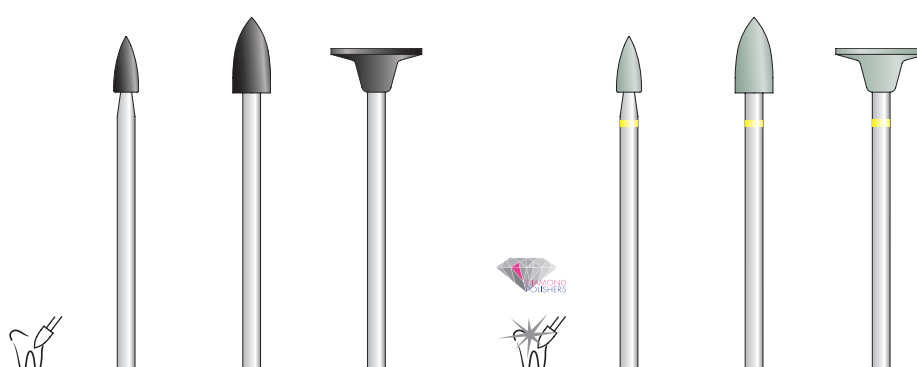


REF	reorder	1006HP	1008HP	1014HP		1106HP	1108HP	1114HP
ISO	652 104 ...	243 524 030	243 524 050	041 524 120	832 104 ...	243 504 030	243 504 050	041 504 120
Size	Ø x L mm	3 x 7	5 x 10	12 x 2		3 x 7	5 x 10	12 x 2
μ	Grit	medium				extra-fine		
	Color	black				green		
rec.		7 500 - 10 000 RPM				5 000 - 7 500 RPM		
max.		10 000 RPM				10 000 RPM		

Box of: 6, 12, 25, 100
N Pressure: 0,3 - 2 N
 Shank: HP, stainless steel

METAL POLISHERS UNMOUNTED


KENDA WHEEL & KNIFE

Gold prepolisher, for gold and semi precious alloys.

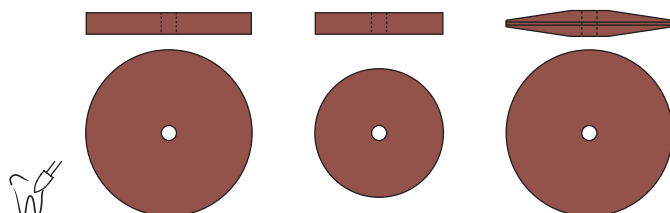


REF reorder	2522-R	2517-R	2522-L
ISO 658 900 ...	372 513 220	372 513 170	303 513 220
Size Ø×L mm	22×3	17×3	22×3
μ Grit	fine		
Color	brown		
rec.	20 000 RPM		
max.	20 000 RPM		

Gold high shine polisher, for hi-shine on gold and semi precious alloys.

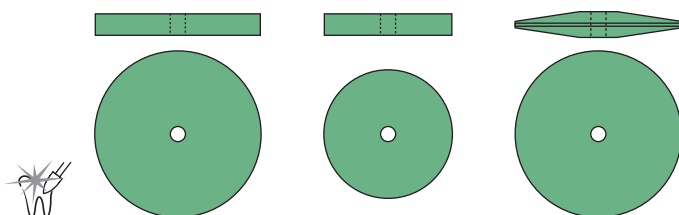

REF reorder	3522-R	3517-R	3522-L
ISO 658 900 ...	372 503 220	372 503 170	303 503 220
Size Ø×L mm	22×3	17×3	22×3
μ Grit	extra-fine		
Color	green		
rec.	10 000 RPM		
max.	20 000 RPM		

Titanium polisher, eliminates scratches and smoothes all titanium surfaces

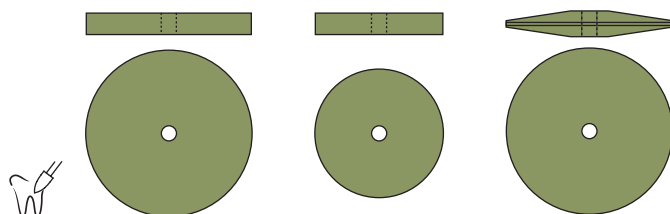

REF reorder	9522-R	9517-R	9522-L
ISO 658 900 ...	372 533 220	372 533 170	303 533 220
Size Ø×L mm	22×3	17×3	22×3
μ Grit	coarse		
Color	oliv		
rec.	20 000 RPM		
max.	20 000 RPM		

Box of: 10, 100
N Pressure: 0,3 - 2 N

COMPOSITE HYBRID POLISHERS

KENDA HYBRID, HP


2-step polishing system

2-step polishing system for all composite materials.

Whitish grey pre-polisher for contouring and smoothing, pale pink diamond polisher for high gloss on composites.



Set REF 791

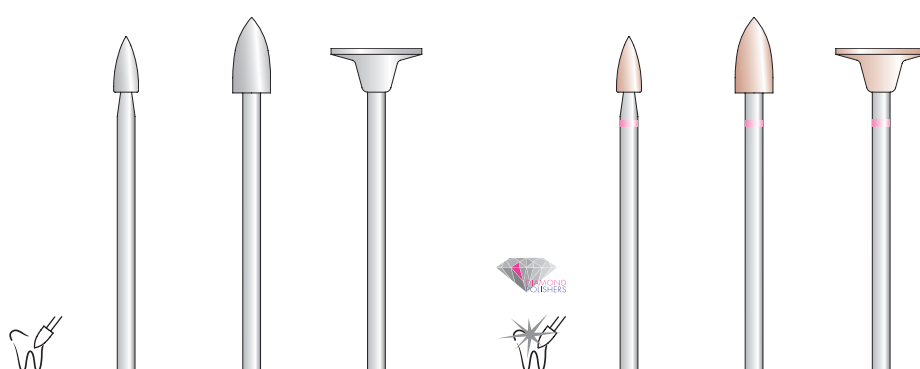

REF	reorder	0006HP	0008HP	0014HP		0106HP	0108HP	0114HP
ISO	632 104 ...	243 524 030	243 524 050	041 524 120	822 104 ...	243 494 030	243 494 050	041 494 120
Size Ø×L mm		3×7	5×10	12×2		3×7	5×10	12×2
μ Grit		medium				ultra-fine		
Color		whitish grey				pale pink		
rec.		7 500 - 10 000 RPM				5 000 - 7 500 RPM		
max.		10 000 RPM				10 000 RPM		

Box of: 6, 12, 25, 100
N Pressure: 0,3 - 2 N
 Shank: HP, stainless steel

PRECIOUS ALLOYS


KENDA UNIVERSAL, HP

2-step polishing system

Ideal for polishing gold, semi precious alloys and composites.
Brown for prepolishing, green for high shine polishing.



	reorder	2106HP	2109HP	2110HP	2111HP	2112HP	2114HP
	652 104 ...	243 524 030	303 524 120	041 524 120	243 524 055	292 524 050	041 524 120
	Size Ø×L mm	3×7	12×3	12×2,5	5,5×13	5×16	12×2
	Grit	medium					
	Color	brown					
	rec.	10 000 - 15 000 RPM					
	max.	15 000 RPM					

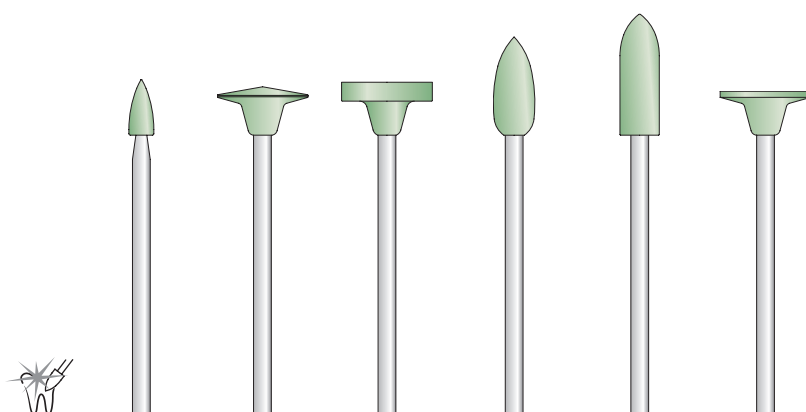

	reorder	3106HP	3109HP	3110HP	3111HP	3112HP	3114HP
	652 104 ...	243 514 030	303 514 120	041 514 120	243 514 055	292 514 050	041 514 120
	Size Ø×L mm	3×7	12×3	12×2,5	5,5×13	5×16	12×2
	Grit	fine					
	Color	green					
	rec.	10 000 RPM					
	max.	15 000 RPM					


Box of: 6, 12, 25, 100
N Pressure: 0,3 - 2 N
Shank: HP, stainless steel

LAB HIGH SHINE POLISHING PASTE

KENDA 244 BLUE

Ideal for high shine polishing of flexible partials as well as full dentures. Also for composites, gold and chrome cobalt, in a stick with approx. 75 ml contents.

For best high shine finish use KENDA SLIMLINE brush made of goat hair bristles and micro fiber leather insole, diameter 48mm, unmounted.

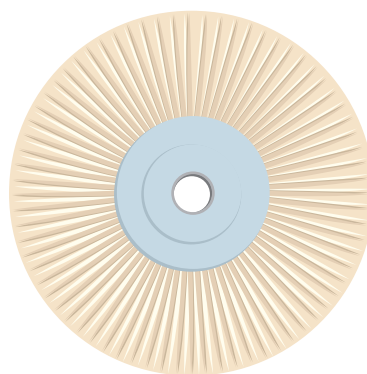

REF reorder	244-BLUE	406
Size	75 ml	48 x 5 mm
Box of:	1, 3	6

AUTOCLAVABLE POLISHER AND BUR CADDY

KENDA PRACTICUS

KENDA PRACTICUS offers space for 6 RA polishers of 22-31 mm length and 12 FG diamond instruments of 16 to 24 mm length.

Thanks to well thought out alignment ergonomic handling with only one hand is possible. Due to the snap lock the caddy can be easily closed and content is fixed for sterilisation. The high-tech medical polymer allows up to 3000 sterilisation cycles in the autoclave at 135°C, 275°F.



REF 270



REF 271



REF 272



REF 273




LABEL

EXPLAINING THE LABEL

The label that you can find on all of our products is designed according to latest standards and has a lot of information on it. It is readable for humans as well as for barcode scanners. The barcode is in the format of «HIBC Data Matrix» which can hold a lot of information. In this example label, the barcode contains our unique HIBC labeler identification code (EKEN), REF, expiration date as well as the LOT number.

Labels can vary for certain products and countries.

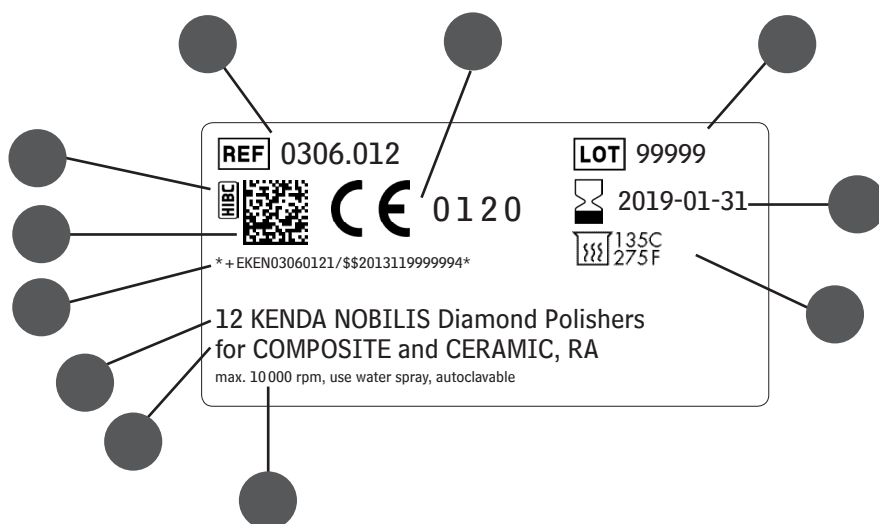

●	order number
●	HIBC logo to identify barcode
●	machine readable barcode
●	human readable barcode information
●	quantity and product name
●	application
●	max. RPM, security advices
●	information for sterilisation
●	expiration date (yyyy-mm-dd)
●	LOT number
●	CE mark with notified body

SYMBOLS

Signs at a glance

Pictographs will appear to indicate the specific use of the instrument.



Polisher for smoothing and finishing



Polisher for finishing and silky shine polishing



High shine polisher for ultimate gloss



autoclavable with moist heat sterilisation



use waterspray for cooling



Order number



LOT number – for traceability of the respective production batch



Expiry date



Manufacturer



user instructions / refer to catalogue



wear eye protection



Diamond particles as material of working part for efficient polishing



maximum speed in RPM. Do not exceed.



recommended speed in RPM for best polishing results.



Disposable. Use one time only.



packaging units in boxes



blister



CE-sign confirms conformity with Directive 93/42/EEC



color



storage temperature



keep away from sunlight



wear breathing protection



KENDA AG • Dental Manufacturing • Industriezone Neugut • Kanalstrasse 18 • LI-9490 Vaduz • Principality of Liechtenstein
Phone +423 388 23 11 • Fax +423 388 23 19 • info@kenda-dental.com • www.kenda-dental.com



© All material herein is copyrighted by KENDA. Printing is permitted without editing. All editing is prohibited.
Contact us for any questions regarding graphical material.