


Helio

Lighting on your hand

Portable Lighting System

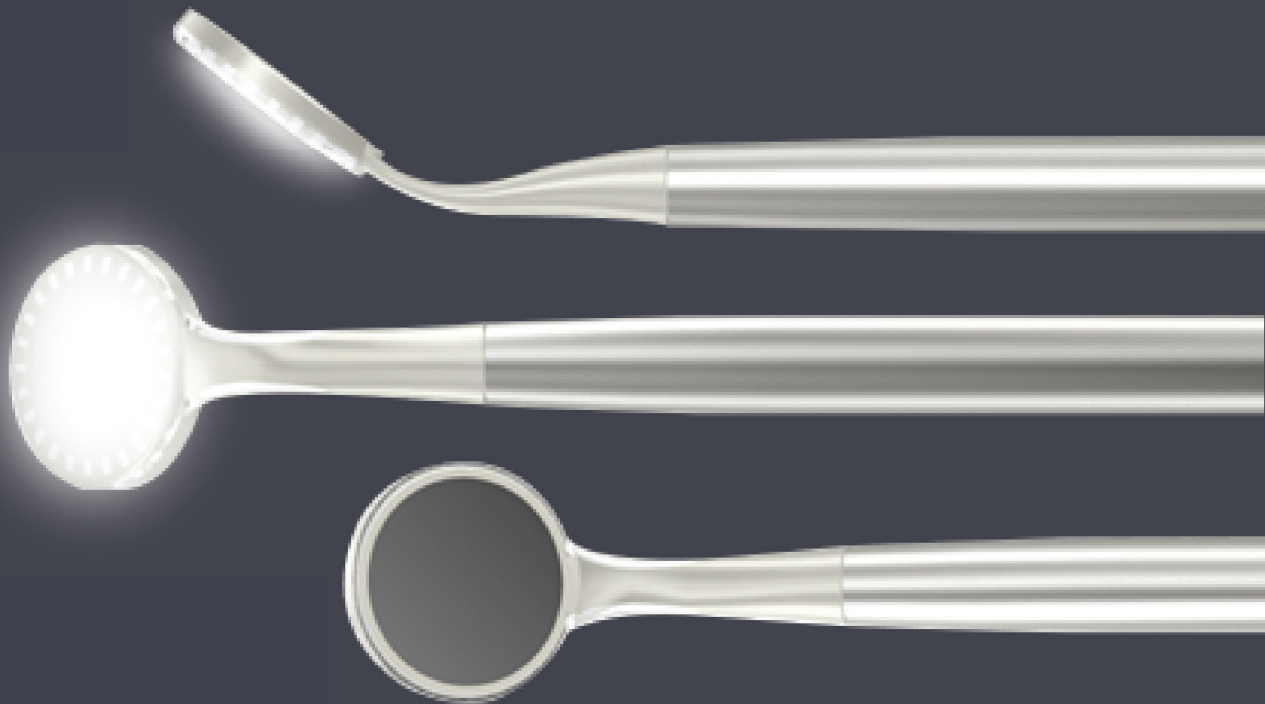

Features

- Lighting anywhere
- Autoclaveable
- Cordless type
- Practical design
- Big capacity battery
- Smart system

Patent pending

Lighting on your hand

Helio is a portable lighting system. Dental mirror shape is most familiar and convenient to make lighting anywhere in oral cavity.

Cordless, Autoclaveable, and many variety functions for users make its value up.

Most Economic, Most Practical Lighting System

There is no switch to turn on/off and change the intensity.

Smart system makes you enjoy using it. Don't need expensive astral lamp, inconvenient head lamp, fiber optic system, or other intraoral lighting system.

Features of Helio



Lighting anywhere

Make light anywhere in oral cavity. Mirror shape is most familiar and easy to use for any kinds of dental practices.



Autoclaveable

The handpiece is autoclaveable. Keep the sterilized condition after autoclave.

⚠ The battery pack is not autoclaveable.



Easy to change Mirror

The mirror interchangeable easily with the magnetic inside of device.



Cordless type

Without wire makes its convenient up.



Practical design

The device can be used like dental mirror.

The device is designed possible to keep sterilized condition for surgical operation.



Big capacity battery

Useable the battery more than 2 hours after full charge.



Smart system

There is no switch. Turn on with grip it, and turn off with put it.

To change intensity, shake the device vertically.

To save battery, the device turn off automatically when don't use it.

Keep & Use Example

Keep

1. Envelop the handpiece in the sterilization pouch.
2. Autoclave the handpiece.

⚠ Do not autoclave the battery pack.

Use

1. Open the sterilization pouch
 2. Insert battery pack
- Be careful to keep sterilization condition.

⚠ The light from device would make polymerization. Please be careful when use resin or other materials.

Example



Specification

| | | |
|--------------|--------|-------------------------|
| Handpiece | Size | Approx. 190 x 25 x 25mm |
| | Weight | Approx. 20g |
| Battery Pack | Size | Approx. 125 x 10 x 10mm |
| | Weight | Approx. 20g |
| Charger | Size | Approx. 105 x 68 x 55mm |
| | Weight | Approx. 100g |

Patent pending